
QUICK START GUIDE

Windows x64 Standard DB Integration on PostgreSQL

Document version 1.00

This guide describes the basics of installing and configuring web-based DICOM Viewer for PACS server. More information is provided in the full Install Manual, download its latest version: [Install Manual](#).

FDA K162011
CLEARED 510(k)




CE 2460
Class IIa certified

ABOUT SOFTWARE


This guide describes MedDream DICOM Viewer deployment on Windows x64. It uses standard database integration on PostgreSQL and can be integrated into any PACS archive, HIS/RIS/EHR/EMR, Telemedicine, Patient Portals, National Systems, CAD and AI algorithms. MedDream Universal Enterprise Viewer can be installed on premises, in any virtual environment or on cloud.

MedDream is designed to aid medical professionals in day-to-day decision-making process. MedDream is a cost effective certified solution perfect for software vendors, integrators, OEM's and national system providers.


PREREQUISITES

-  The Operating system is **Windows 10 (x64)** or one of the **Windows Server (x64)** versions: 2012, 2012 R2, 2016, 2019.
-  **PostgreSQL database** is installed. If PostgreSQL is not installed on the server, you can download it here: <https://get.enterprisedb.com/postgresql/postgresql-10.8-4-windows-x64.exe>
During the PostgreSQL server configuration remember the port and superuser password.
-  Make sure the server has at least **2 GHz CPU with 2 cores, 3 GB of RAM, 20 GB** of available **hard disk** space. Network bandwidth is not less than **4 Mbit/s**.


INSTALLATION FOR THE FIRST TIME

-  **Install MedDream DICOM Viewer:**
 - **Download** the latest Windows version of MedDream from the Softneta website: <https://www.softneta.com/files/meddreamviewer/760/MDPACK-7.6.0.exe>
 - **Execute MDPACK-7.6.0.exe.**


- Remember default MedDream Viewer **login credentials** and click **Next >** button:

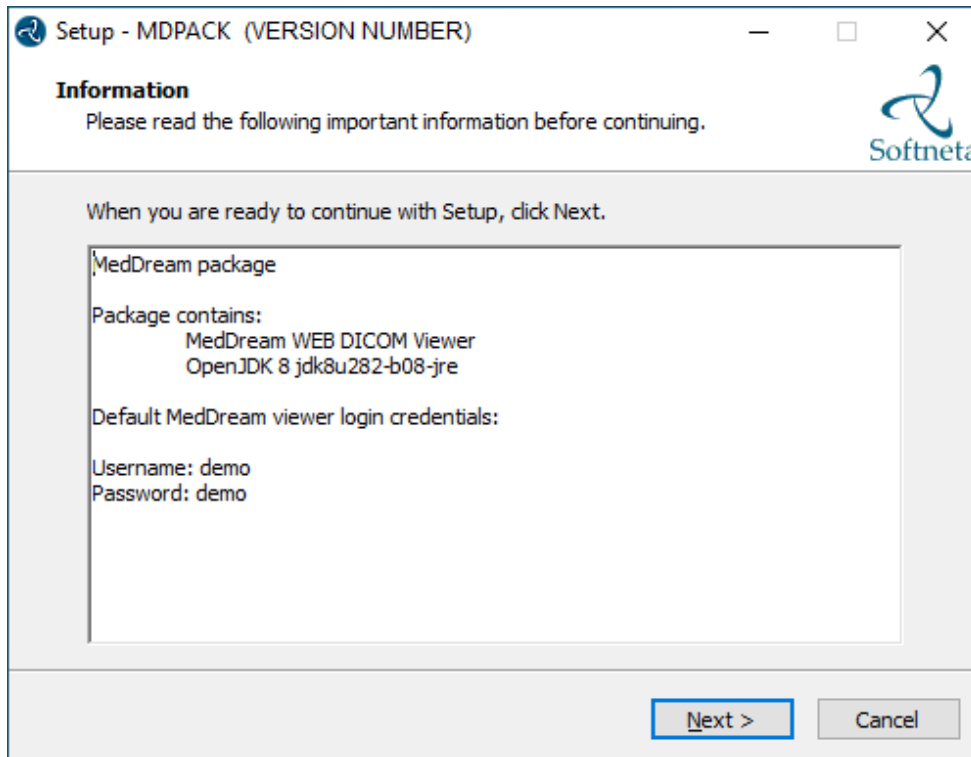


- Choose **C:\MDPACK** and click **Next >** button:

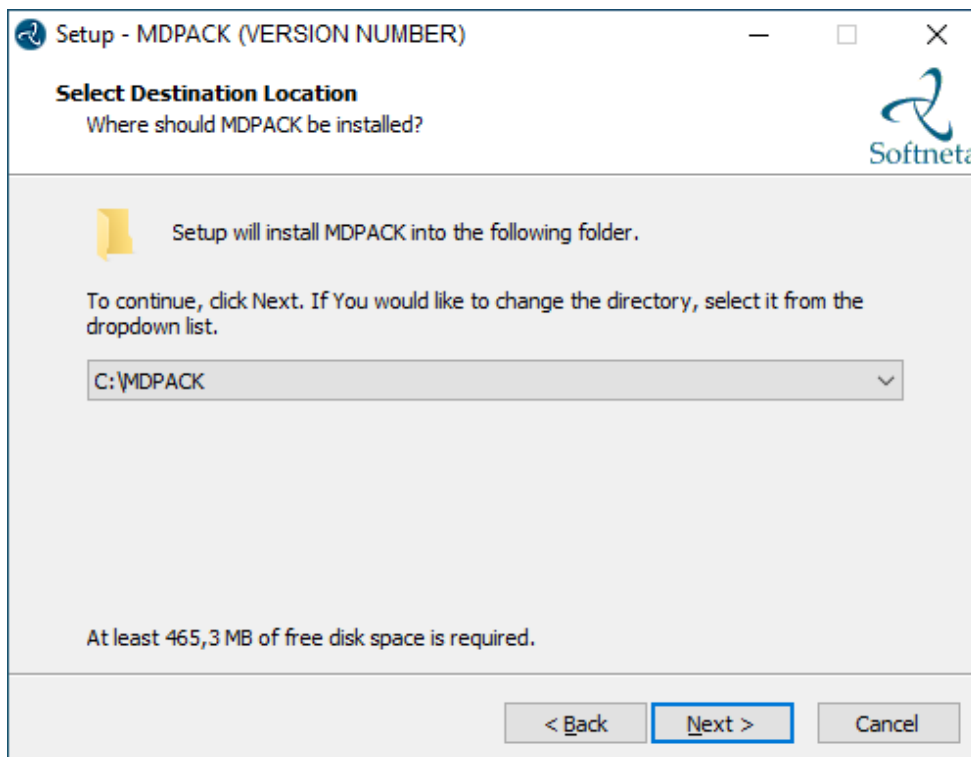


- Choose **Default installation** and click **Next >** button:

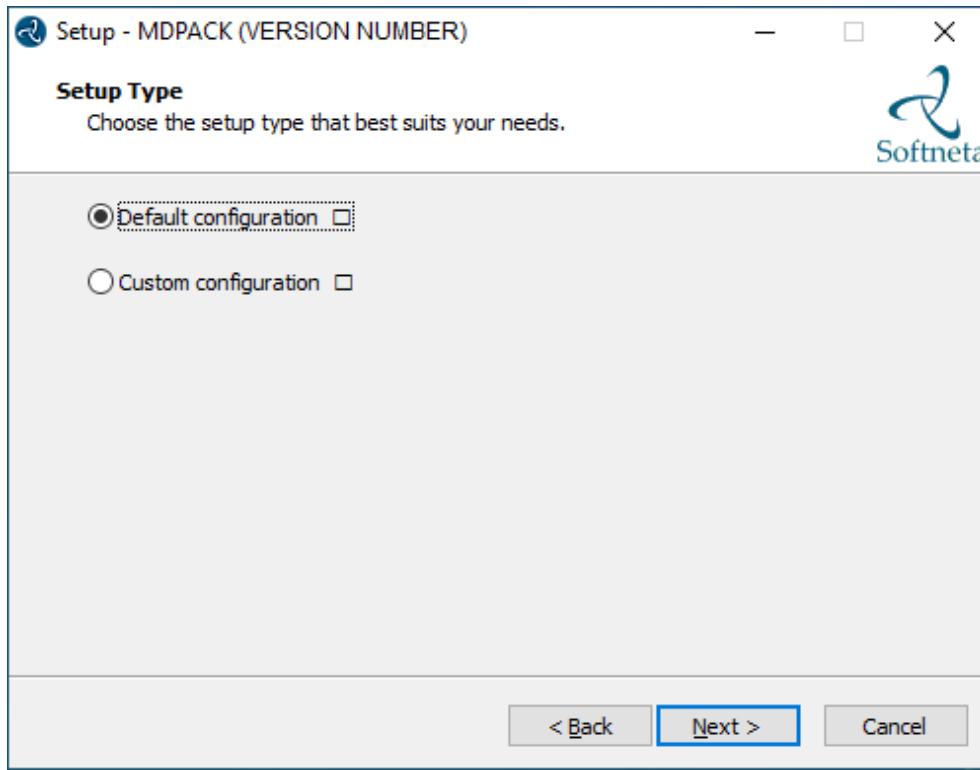

- Click **Install** button:

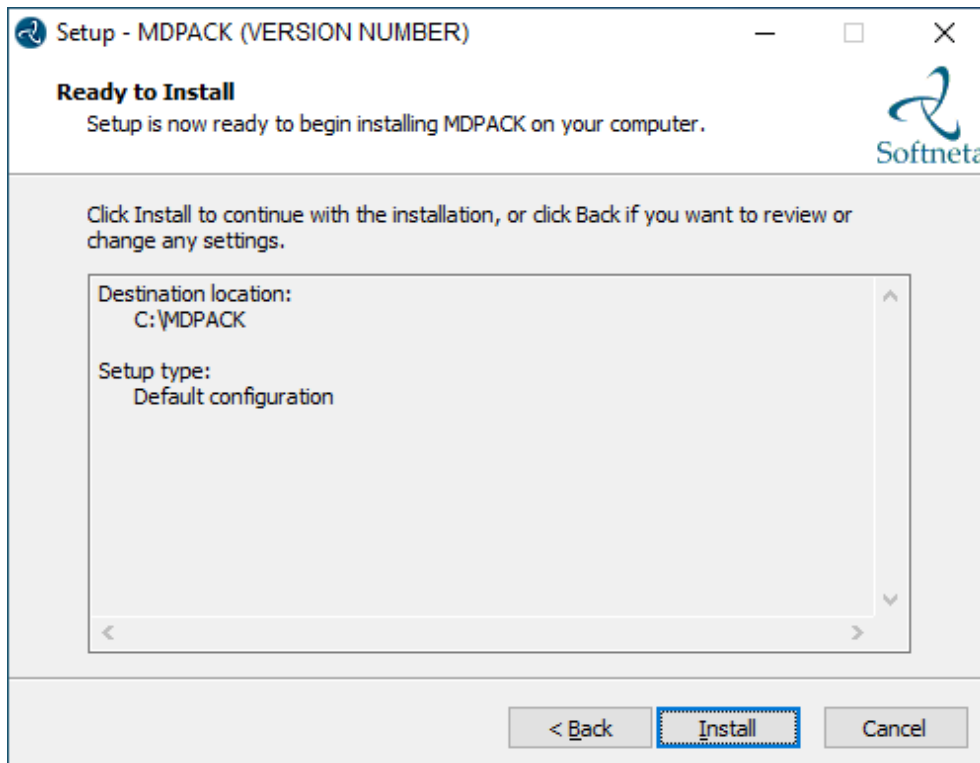


- Well done, you have completed the installation. Click **Finish** button.

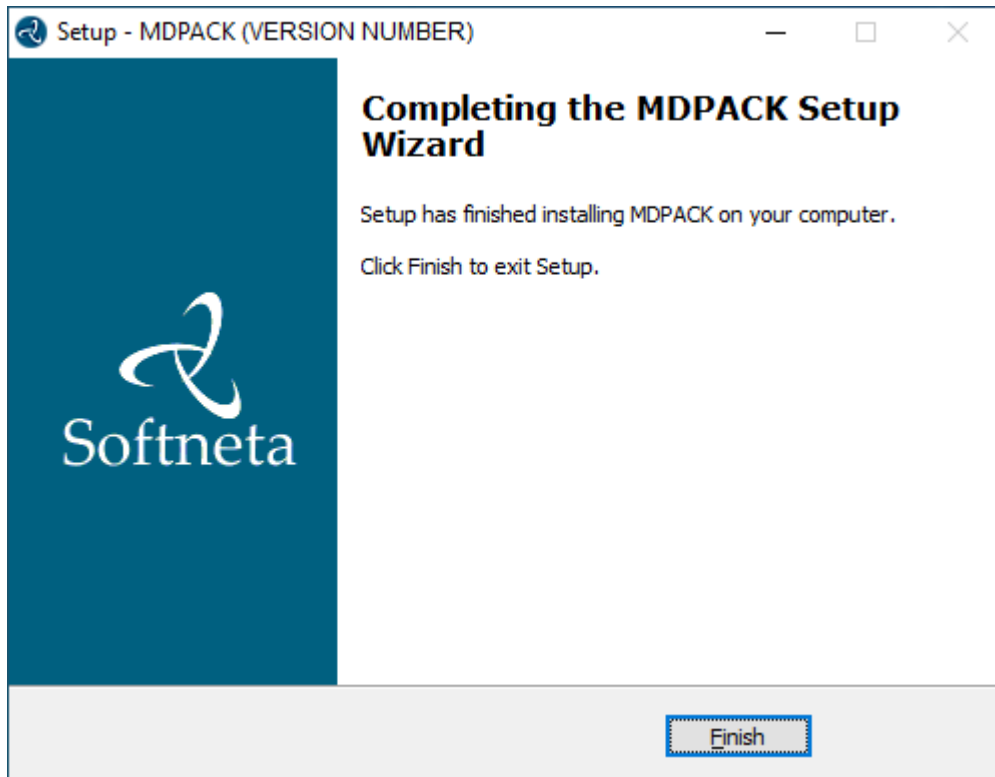


2 Install MedDream database:

- **Download** scripts: <https://www.softneta.com/files/meddreamviewer/752/MedDreamViewer-7.5.2-DB-PostgreSQL.zip>
- **Unzip** the package MedDreamViewer-7.5.2-DB-PostgreSQL.zip to **C:\MDPACK\MedDreamViewer-7.5.2-DB-PostgreSQL**
- **Install** MedDream database to PostgreSQL by executing C:\MDPACK\MedDreamViewer-7.5.2-DB-PostgreSQL\PostgreSQL-DB\meddream-db-create.bat:

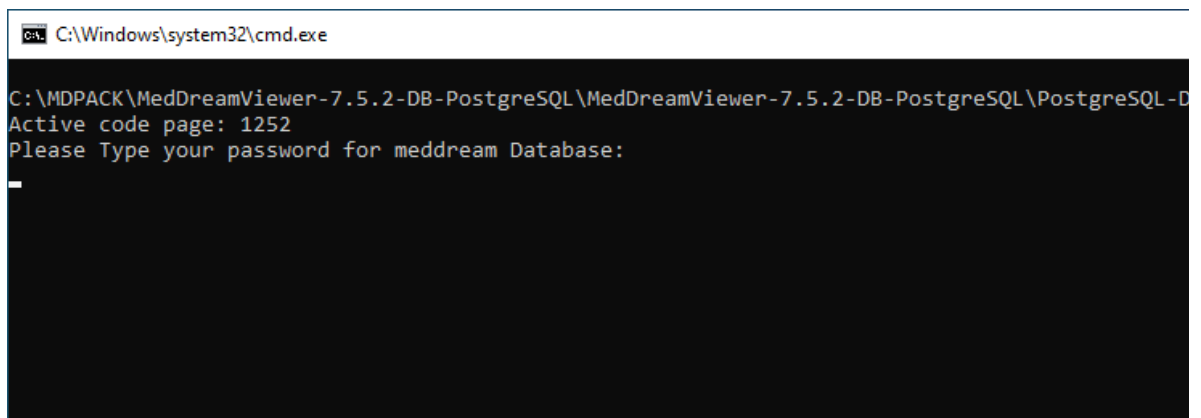

- **Type** your password and press **<Enter>**. This is a new password for the user "meddream". **Remember it.**
After that you will be asked to enter the current password for the user "postgres". **Type** "postgres" user password of **PostgreSQL database** and press **<Enter>**.
- Insert test data to MedDream database by executing C:\MDPACK\MedDreamViewer-7.5.2-DB-PostgreSQL\PostgreSQL-DB\meddream-db-insert.bat:

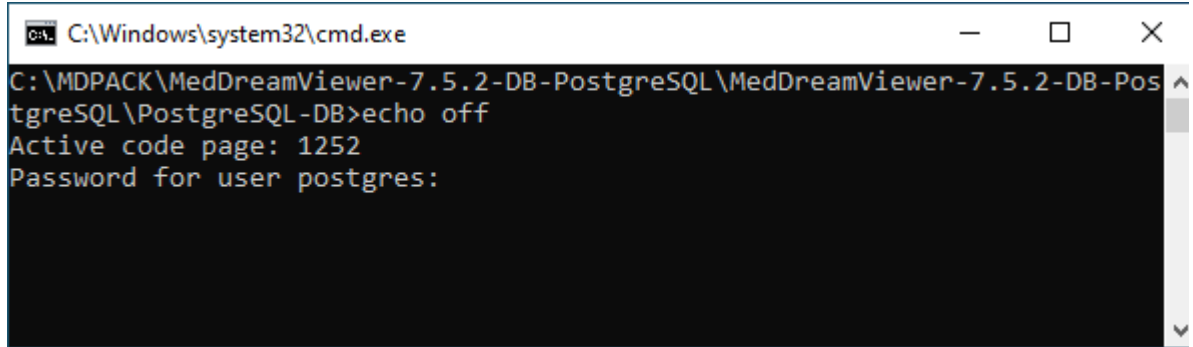

```
C:\Windows\system32\cmd.exe
C:\MDPACK\MedDreamViewer-7.5.2-DB-PostgreSQL\MedDreamViewer-7.5.2-DB-PostgreSQL\PostgreSQL-DB>echo off
Active code page: 1252
Password for user postgres:
```

- You will be asked to enter the password for the user "postgres". **Type** it and press **<Enter>**.

CONFIGURATION

1

Copy and replace C:\MDPACK\MedDreamViewer-7.5.2-DB-PostgreSQL\MDPACK folder to C:\MDPACK

2


Change the basic MedDream Viewer settings for the first launch.

- **Open** C:\MDPACK\MedDream\application.properties file and **write** database user "meddream" password after "=" (you created this password during MedDream database installation):
com.softneta.meddream.pacs.configurations[0].password={Password} #line 418
- **Save** the file.


3

Restart MedDream service. After making any changes to the application.properties file, the MedDream service must be restarted.

- **Open** Windows **Services** by pressing **<Win> + R** and **type** services.msc. Click **<OK>** button:

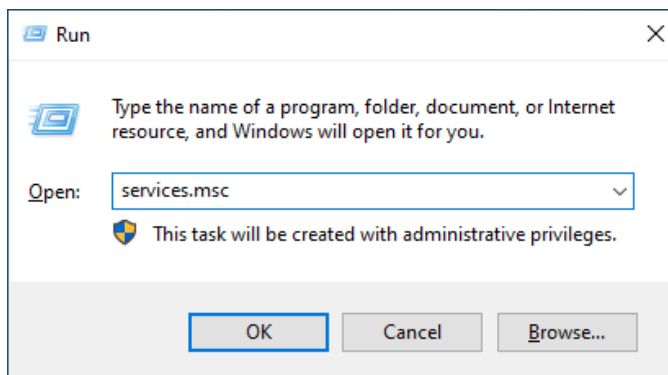

- Find the **MedDream** service and restart it:

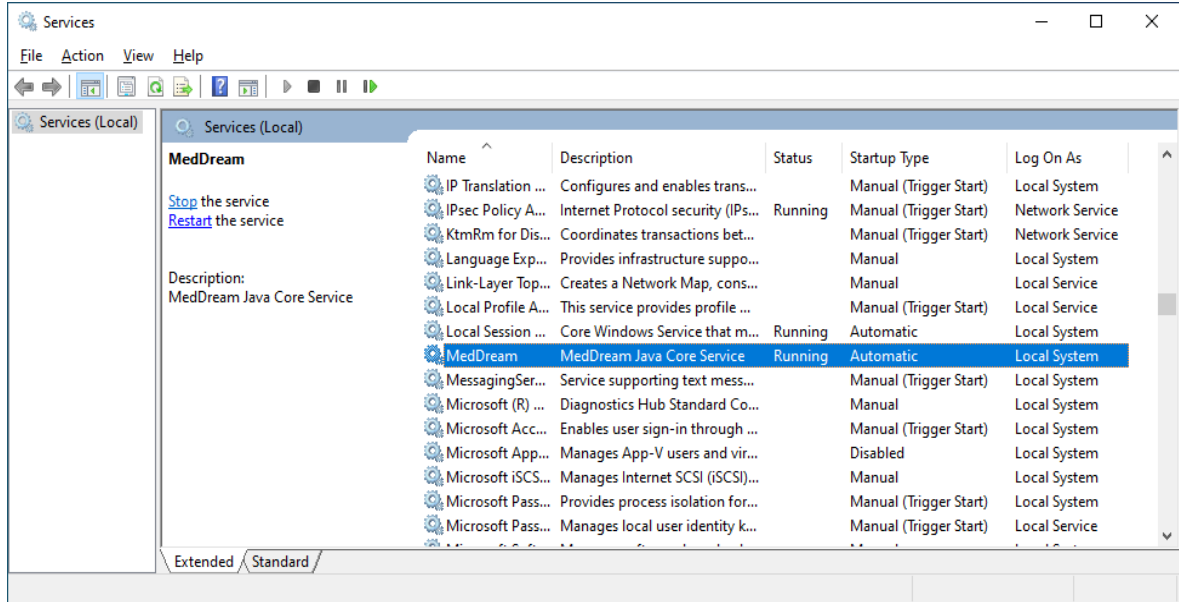

NOTE: In the C:\MDPACK\MedDream\ApTemplates you can find examples of MedDream Viewer configurations with explanations. All MedDream settings are stored in the application.SAMPLE.properties file.


VERIFICATION

1

Open a browser and enter the MedDream URL <http://127.0.0.1/> in the address bar.

2


In the login window, enter the username **demo** and password **demo**. Press **Log in** button:

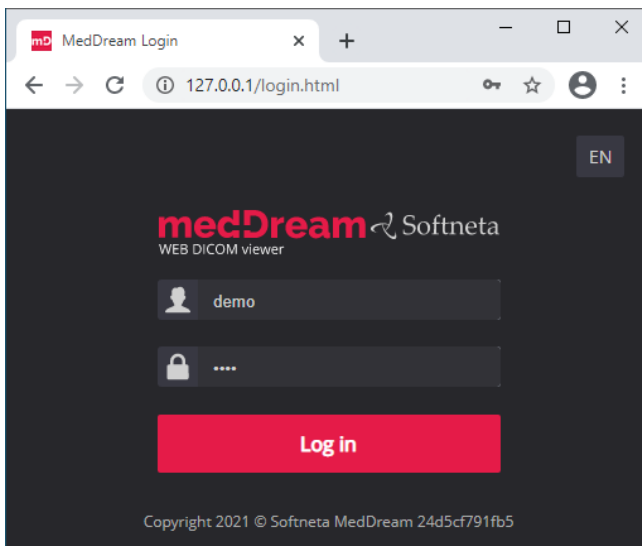

TIP: If you can't connect, restart Windows services (run in cmd line "services.msc" and restart MedDream service).


3

MedDream Viewer shows an empty list of studies:

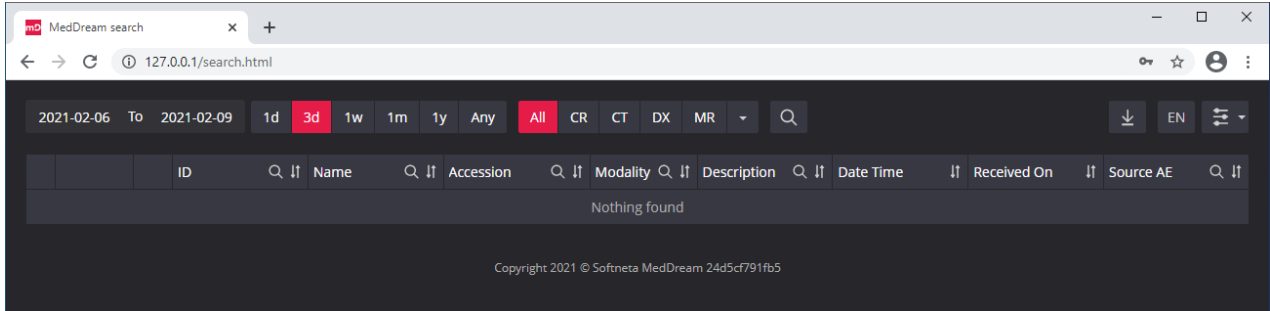

4


Change the date range selection by pressing **<Any>** button. Now you can see one study in the list:

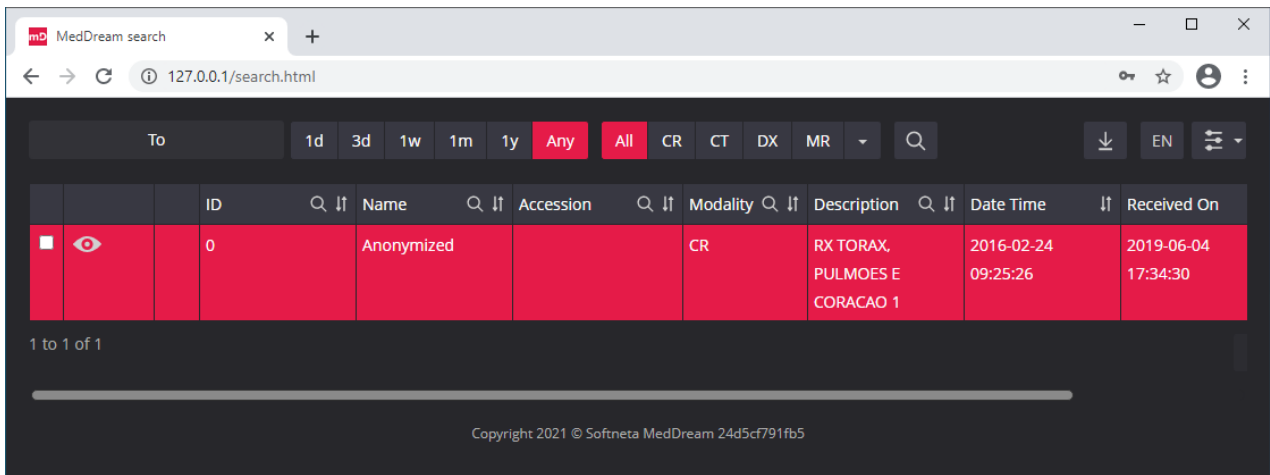

TIP: Samples of anonymized studies can be downloaded from <https://demo.softneta.com/> (Login → Open the selected study → Click Export → Select the DICOM format and click Export again → Press the Download link).

5

Click once on the study with the left mouse button. MedDream viewer opens the study in a new browser tab:

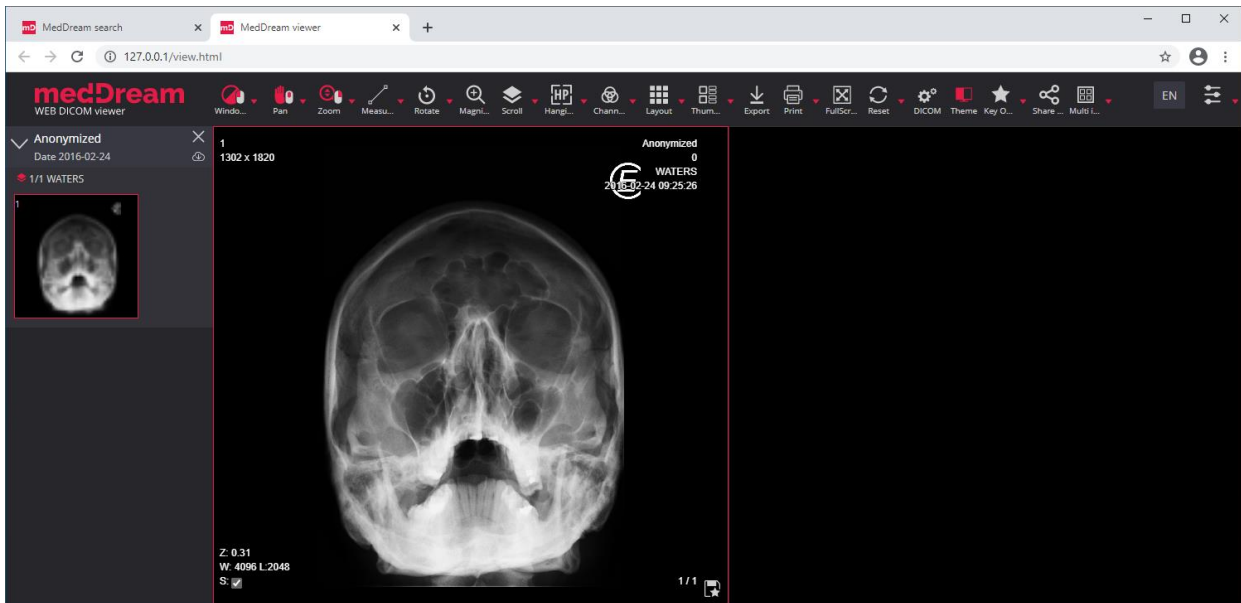

6

Now you can try and see how powerful and useful this tool is.

TROUBLESHOOTING

✓

You must restart the MedDream service each time after changes to the application.properties file or if something went wrong. You can restart the service by running in cmd line: `services.msc` -> then selecting MedDream service and clicking Restart on the left.

✓

If you can't open the studies, clearing your browser's cache may help. Open a browser and press **CTRL+SHIFT+DEL**. Then select which old data you want to delete and click **Clear** data.

✓


MedDream DICOM Viewer logs are stored in the **C:\MDPACK\MedDream\meddream.main.log**. The standard logging level in the application.properties file is set to INFO. If the log contains too little information, you can change the logging level to **DEBUG** by changing the value of the logging.level.com.softneta parameter.

LICENSE

MedDream DICOM Viewer licensed under the Commercial License. Please contact info@softneta.com to get TRIAL or Commercial license. If you have a License Serial number, follow these instructions:

1

Press **Register** button. The system opens a license registration window:

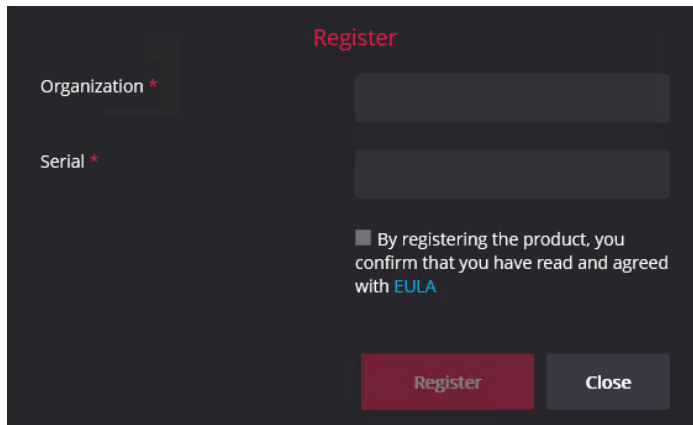

2

Enter the name of the organization and the license **Serial** number you received.

3

Click on the link **EULA** and read Software License Agreement. After reading, close the window by clicking the **<Close>** button.

4

Check the box to confirm that you accept license agreement.

5

Press **Register** button. The license will be successfully registered.

NOTE: If you have any questions about MedDream licenses and prices, please contact our SALES team by e-mail info@softneta.com or by phone +370 670 15991.

HELP AND RESOURCES



Detailed MedDream installation and configuration instructions are provided in the install manual. You can download it from here: <https://www.softneta.com/files/meddreamviewer/doc/MedDream-DICOM-Viewer-Install-Manual.pdf>



We are constantly improving this product. Get the most out of your MedDream Viewer with updated features and performance.



Download the User manual from <https://www.softneta.com/files/meddreamviewer/doc/MedDream-DICOM-Viewer-User-Manual.pdf> and get more information about this device.



Read short instructions and watch tutorial videos here: <https://www.softneta.com/documentation/user-guide/user-introduction/>



Our great SUPPORT team is always ready to help you. If you have any questions, please write an email to support@softneta.com or call +370 645 06960.