


# SendToPACS

CONVERT NON-DICOM  
FILES TO DICOM  
AND SEND TO PACS



VERY EASY  
TO USE

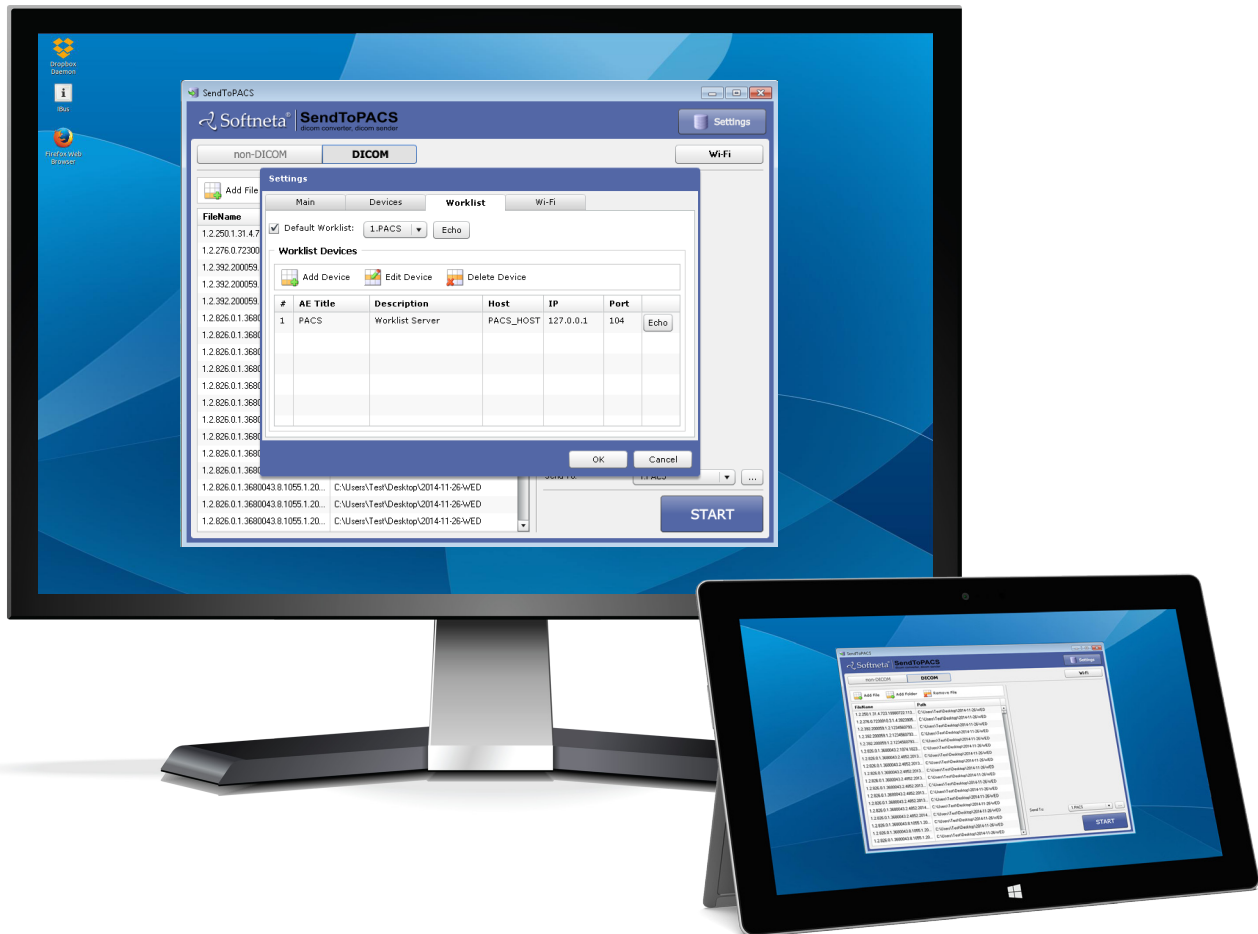

IMAGES, VIDEOS,  
SIGNALS AND  
DOCUMENTS  
SUPPORT



ENHANCED  
CHARACTER SET  
SUPPORT




DICOM MODALITY  
WORKLIST



 Softneta

## Non-DICOM files to PACS

SendToPACS software converts non-DICOM images and video files into DICOM format files, and sends them to PACS server or DICOM workstation.

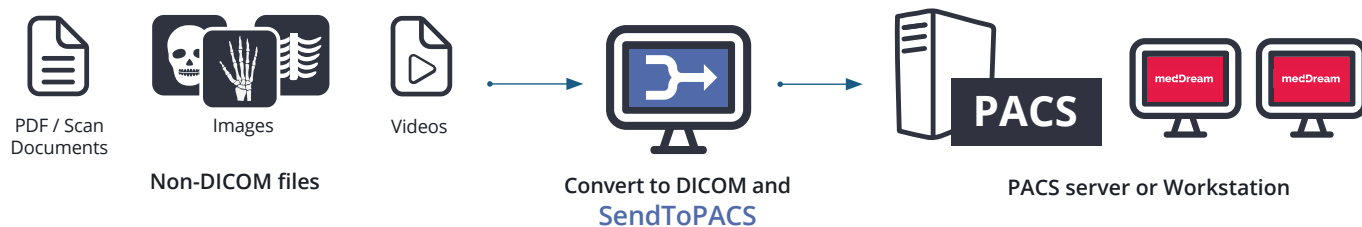

## DICOM files to PACS

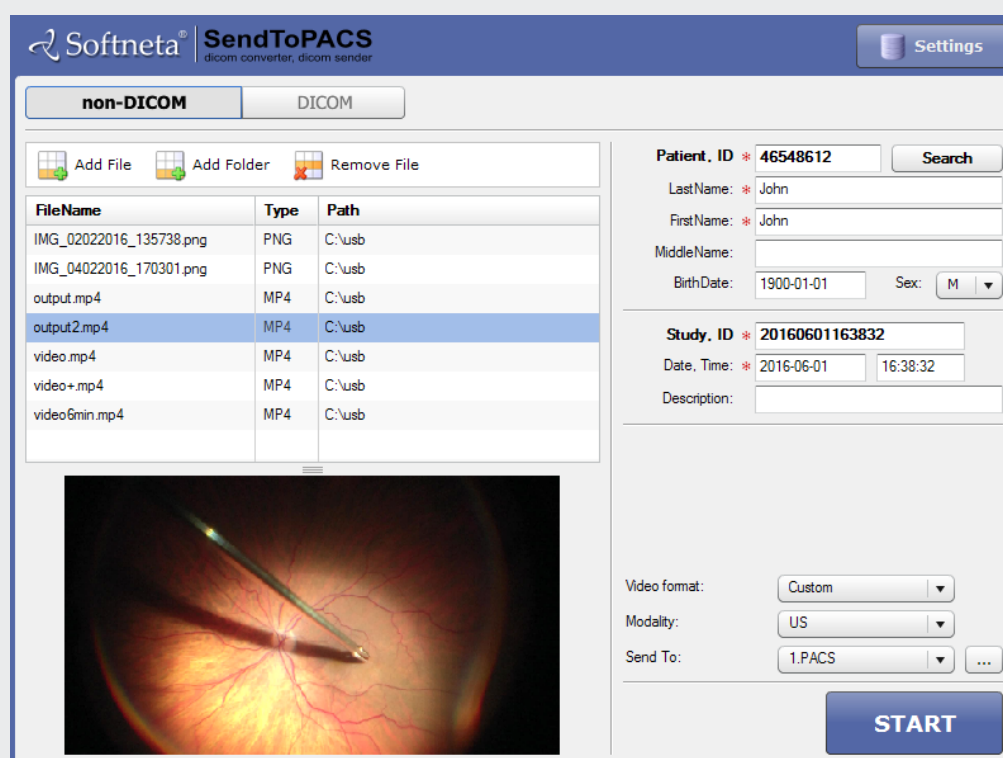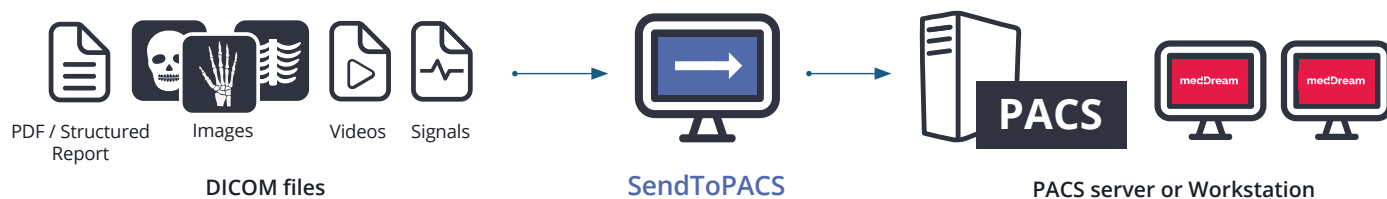

## SendToPACS features:

### Main features:

- Files conversion to DICOM with the highest possible quality;
- Selectable video conversion to DICOM as MPEG2 or MPEG4 format;
- Sending DICOM files to PACS server or DICOM workstation;
- DICOM modality worklist support, advanced worklist search;
- Connectivity verification (echo);
- Windows context menu support;
- Drag and drop files or folders;
- Multi-language support.

### Customization:

- Displays worklist result in UTF-8 or Latin1;
- Possibility to query HTTP/HTTPS worklist;
- Customizable modality representation value;
- DICOM specific character set latin1 (ISO\_IR 100) or UTF-8 (ISO\_IR 192);
- Customizable SOP Class UID;
- Customizable video conversion process CPU resources (number of threads).

### Installation requirements:

- Supports 32-bit & 64-bit Windows OS.

### SendToPACS supports the following files and formats

#### Images:

jpg, png, tif, bmp, dcm.

#### Videos:



mov, mpeg, mpg, avi, mp4, mov, wmv.

#### Documents:

Structured Report and PDF files.

#### Signals:

ECG and others.

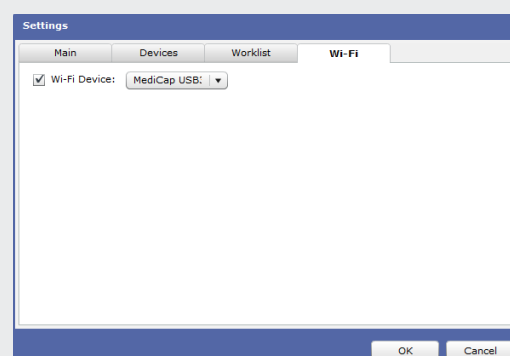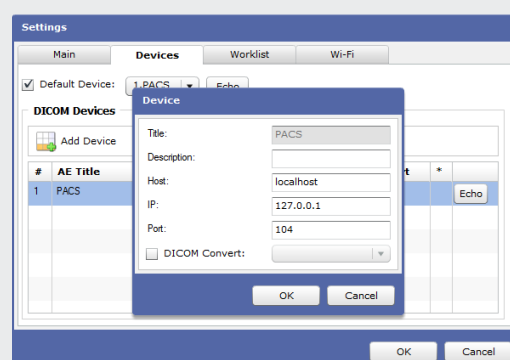

## Workflow with SendToPACS

